


Your
IVD Specialists!

DIA.
PRO

CLIA

A **chemiluminescence** technology
based on **magnetic particles**

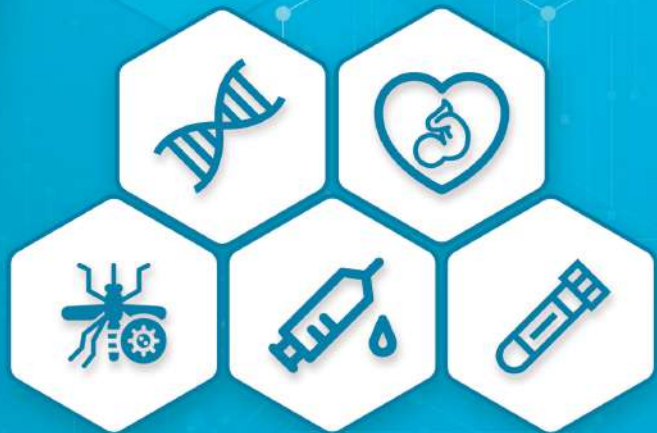

✉ INFO@DIAPRO.IT
🌐 WWW.DIAPRO.IT
📍 VIA GIOSUÈ CARDUCCI, 27, 20099, MILAN (MI) - ITALY



A RAPIDLY EXPANDING MENU



BLOOD SCREENING MARKERS

HEPATITIS B VIRUS

HBs Ag - [RASAG.CE](#)
HBc Ab - [RABCAB.CE](#)

HEPATITIS C VIRUS

HCV Ab - [RACVAB.CE](#)

HIV

HIV Ab&Ag - [RAIVCOMB.CE](#)

HTLV

HTLV I & II Ab - [RAHTLVAB.CE](#)

SYPHILIS (TREPONEMA PALLIDUM)

Syphilis Ab - [RASIAB.CE](#)

CHAGAS (TRYPANOSOMA CRUZI)

T.Cruzi Ab - [RATCAB.CE](#)



EMERGING TROPICAL DISEASES

DENGUE VIRUS

Dengue virus Ab - [RADENGVAB.CE](#)

CHIKUNGUNYA

Chiukungunya Ab - [RACHIKVAB.CE](#)

ZIKA VIRUS

Zika Ab - [RAZIKVAB.CE](#)



SEROLOGY

TETANUS

Tetanus toxoid IgG - [RATETG.CE](#)

HEPATITIS E VIRUS

HEV Ab - [RAEVAB.CE](#)

HEPATITIS D VIRUS

HDV Ab - [RADAB.CE](#)

MALARIA

Malaria Ab - [RAMALAB.CE](#)



PRENATAL SCREENING

TOXOPLASMA GONDII

Toxoplasma IgM - [RATOXOM.CE](#)
Toxoplasma IgG - [RATOXOG.CE](#)

RUBELLA VIRUS

Rubella IgM - [RARUBM.CE](#)
Rubella IgG - [RARUBG.CE](#)

CYTOMEGALOVIRUS

Cytomegalovirus IgM - [RACMVM.CE](#)
Cytomegalovirus IgG - [RACMVG.CE](#)

HERPES SIMPLEX VIRUS

HSV1 IgM - [RAHSV1M.CE](#)
HSV1 IgG - [RAHSV1G.CE](#)
HSV2 IgM - [RAHSV2M.CE](#)
HSV2 IgG - [RAHSV2G.CE](#)
HSV1&2 IgM - [RAHSV1M.CE](#)
HSV1&2 IgG - [RAHSV1G.CE](#)



AUTOIMMUNITY

IgG anti-SSB - [RASSB.CE](#)
IgG anti-SSA52 - [RASSA52.CE](#)
IgG anti-RNP - [RARNP.CE](#)
IgG anti-Sm - [RASM.CE](#)
IgG anti-JO1 - [RAJ01.CE](#)
IgG anti-SCL70 - [RASCL70.CE](#)

Kit format: 100 tests



JARA

A CHEMILUMINESCENCE SYSTEM
DESIGNED FOR EXCELLENCE



Sample Disposable Tips
(Cost effective)
& disposable cuvettes

**Refrigerated Reagent
loading & Storage Area**

2D Matrix Code label for:

- Assay identification
- Lot identification
- Expiration date
- Calibration values



User friendly Touch Screen

**Barcode Scanner Sample
& Control Identification**



Samples & Controls Management

Continuous loading
90 samples and control slots
Sample volume **10-200 µl**
Primary tubes (100-160 mm)

No carry over
Clot & sample volume detection
Barcode **sample** and **control**
identification



ASSAY CARTRIDGE:

- The two cartridges included in the kit contain reagents for 2x50 tests
- Include calibrators and controls
- All reagent vials are sealed with caps to avoid spilling
- Stability of 4 weeks once opened



JARA



**FULLY AUTOMATED BENCH TOP ANALYZER
WITH CONTINUOUS RANDOM ACCESS**

TECHNICAL SPECIFICATIONS

Throughput	Up to 120 tests/hour. First result after about 15 mins. (assay dependent) Then 1 result/in 30 secs.
Mode of operation	Random Access including Continuous Loading and STAT
Sample & Control identification	Barcode reader for samples and controls
Sample dispensation	Standard universal disposable tips to avoid carry-over Liquid level sensor
Reagent	N° 16 cartridges on-board simultaneously
Calibration	Calibrator included into the cartridge for assay calibration Calibration Curve memorized on RFI0/20 barcode label
Controls	Ready-to-use liquid, included in the kit box
Operating System	Touchscreen on board Password for system access Color-coded graphical user interface Ergonomic flat, low-glare, LCD, touchscreen
Configuration	Bench top, compact design
Dimension and weight	Heigh 90 cm – Width 190 cm – Depth 78 cm – Weight +/- 200 kg
Certification	CE (IVD Directive 98/79)



SARA

S.A.R.A. CONSUMABLES

CODE	ITEM	DESCRIPTION
RALPACK.CE	<i>DIA.CHEMILUX Luminol Pack</i>	Kit of 3 x 150 mL Lum.A and 3 x 150 mL Lum.B
RAWSPACK.CE	<i>DIA.CHEMILUX Washing Pack</i>	20x Concentrated (30 x 50mL)
06311	<i>SARA Starter Kit</i>	Kit for setup and calibration of SARA Instrument
03575	<i>SARA Tips Pack</i>	Minimum Pack Format - 1 Box -20 x 500 pcs (10 000 pcs)
06352	<i>SARA Cuvette Pack</i>	Minimum Pack Format - 90 x 98 pcs/box (8820 pcs)
06218	<i>SARA System Liquid</i>	Needles Solution Pack for SARA Instrument
06203	<i>SARA Essential Spare Part</i>	Minimal required Spare parts
06219	<i>SARA Decontamination Kit</i>	Kit for decontamination of SARA Instrument
06502	<i>SARA Maintenance kit</i>	Kit for maintenance of SARA Instrument



For more information visit:
WWW.DIAPRO.IT



Dia.Pro Diagnostic Bioprobes S.r.l

Via G. Carducci, 27 - 20099
Sesto San Giovanni (MI) - Italy
Tel +39 02 27.00.71.61
www.diapro.it
email: info@diapro.it

EN ISO 13485 Certified
BRO-CLIA Rev0923

